
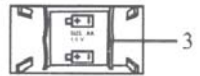
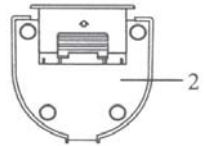


# INSTRUCTIONS

Please read the instructions carefully before use.

1. Denture Cleaner requires 2 x "AA size" batteries, 1.5V, (not included).
2. Take out the cleaning compartment (fig. 1) and wash with clean water.
3. Wipe the cleaning compartment with a cloth and then place it back on the base unit.
4. Open the battery compartment (fig. 3) at the bottom of the base (fig. 2) and insert 2 x "AA size" batteries according to battery + / - direction. Close the lid properly.
5. Fill the cleaning compartment with water or your own denture cleaner.
6. Place the denture(s) into the cleaning compartment, then turn the machine on for 15 minutes. (If denture(s) is stained, light brushing may be necessary).
7. 5 minutes of daily denture cleaning is recommended. If your denture(s) has accumulated stain and tar deposits, then clean for 15 minutes daily until stain and tar are removed.
8. After each cleaning, use warm water to rinse your denture(s) thoroughly.
9. Wash the cleaning compartment with water (fig. 1) and dry with cloth after use.
10. TURN OFF the denture cleaner after use.

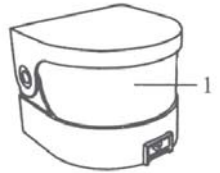

## NOTES

- Before placing denture(s) in the cleaning unit, remove any adhesive.
- This unit will work with any liquid denture cleaner.
- Please do not immerse the main unit in the water.
- Please replace with new batteries if the power gets low.
- Do not mix different battery types together.
- When replacing batteries, replace all batteries at the same time.
- When not in use for a long period, take out the batteries to avoid battery oxides which will damage the unit and impair its function.