

Mimitakara[®]



@goodHearing A1

User Guide


1. ACCESSORIES	P 1
2. FEATURES	P 2
3. SET UP	P 3
3.1) Wearing Instructions	P 3
3.2) Recharging Instruction	P 5
3.3) Switching Between Modes	P 7
4. PAIRING BLUETOOTH	P 9
4.1) Download App.....	P 9
4.2) Bluetooth Connect	P 9
4.3) Notes	P 12
5. SAFETY & MAINTENANCE	P 14
5.1) Earbud Maintenance	P 14
5.2) Battery Maintenance	P 15
5.3) LED Statuses	P 16
6. FAQ	P 18
7. SPECIFICATION	P 21
8. WARNING AND REMINDER	P 22

1. ACCESSORIES

- ① Main device with replaceable earbuds (M size)
- ② 3 x Additional pair of L, S & Mini size earbuds
- ③ Cleaning Brush
- ④ USB Cable



2. FEATURES

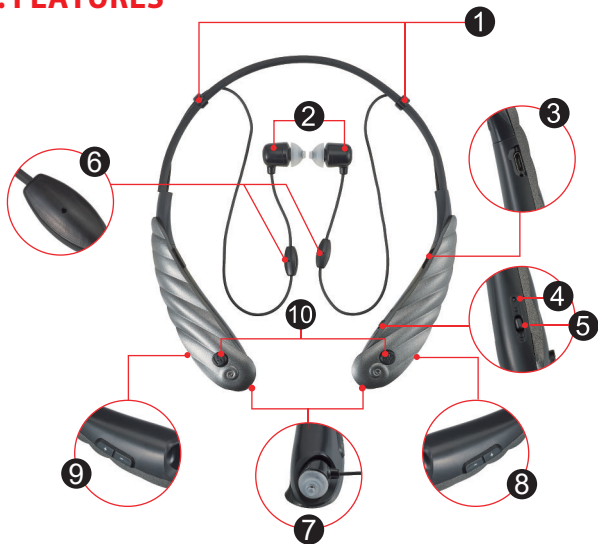

- ① Cable Clips
- ② Earphones
- ③ Charging Port
- ④ LED Indicator
- ⑤ Power Switch
- ⑥ Dual Microphone
- ⑦ Magnetic Earphone Holders
- ⑧ Left Volume Controller
- ⑨ Right Volume Controller
- ⑩ Mode / Bluetooth Connect

3. SET UP

3.1) Wearing Instructions

- 1) Simply place the device around your neck, under your shirt collar or on the neck and place the earbuds into your ear.



- 2) The power switch is on the left side of the device facing your neck.



3) The **goodHearing A1** has adjustable volume settings. To adjust the volume of the device, the right and left sides **+** / **-** buttons situated on the outer sides base of the device, increase and decrease the volume respectively to suit your preferred specification

CAUTION: LISTENING AT A HIGH VOLUMES MAY DAMAGE YOUR HEARING

4) We suggest placing the device around your neck, under your collar to avoid it from rubbing against objects in order to reduce static sounds (Ensure the microphone is not obstructed for optimum sound quality)

3.2) Recharging Instruction

- 1) Charging the battery before use is recommended but not required. The **goodHearing A1** comes pre-charged however initial charge will vary.



- 2) 3+ hours for full charge gives up to 10 hours of usage time



3) When plugged into a power source



Solid RED LED Indicator : Charging



Solid BLUE LED Indicator : Fully charged

3.3) Switching Between Modes

1) There are 3 pre-set modes on the **goodHearing A1**.

- General : Standard Mode


Mostly used indoors where the environment is not too noisy.

- Conversation : Noise Cancelling Mode

Used in louder environment such as outdoor public areas, restaurants.

- In-car : High-level Noise Cancelling Mode

Providing superior noise cancelling in even the loudest environments, such as in transportation.

2) To switch between the 3 modes click the  button (located on the center of the **goodHearing A1** device). A beep will be heard when switching between the 3 modes to indicate which mode it is currently set.

- Mode 1 : A single beep will be heard
(General: Standard Mode)
- Mode 2 : Two beeps will be heard
(Conversation: Noise Cancelling Mode)
- Mode 3 : Three beeps will be heard
(In-car: High-level Noise Cancelling Mode)

4. PAIRING BLUETOOTH

4.1) Download App

Please scan QR code to download the App or search “goodHearing” on App Store / Google Play



goodHearing




iOS




Android

4.2) Bluetooth Connect


1. Turn on **goodHearing A1** and have both devices within arm length of each other.
2. Press and hold the **(M)** for 3 seconds, the LED Indicator will flash alternating BLUE and RED colors.
3. Turn on the Bluetooth function on your mobile device.

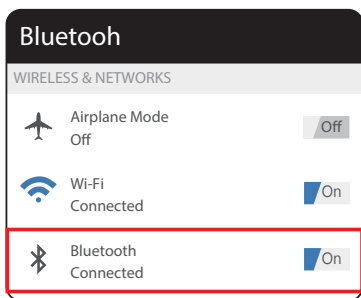

4. Scan for nearby devices, search and pair with **goodHearing A1**.

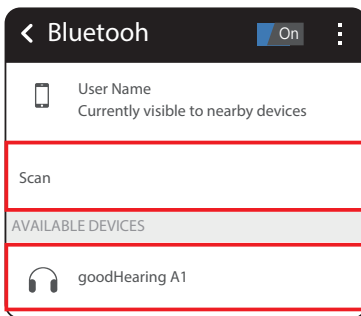

5. When pairing has been successful your phone will show that the status is “Connected” with **goodHearing A1**.

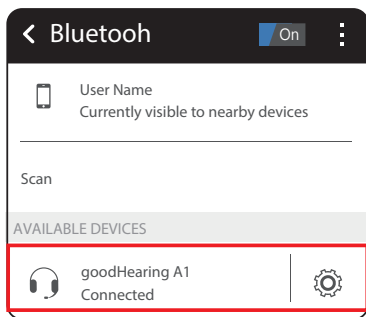

- ✂ The pairing steps may vary according to your mobile operating system.

4.3) Notes

- Paired devices will remember each other unless pairing information is deleted manually, so there is no need to pair again. The **goodHearing A1** can pair and remember 2 Bluetooth devices.
- Repeat the steps outlined above to set-up your **goodHearing A1** with another Bluetooth device.
- When a 3rd Bluetooth device is paired, the oldest pairing device information with **goodHearing A1** will be deleted.
- The **goodHearing A1** will connect automatically with the latest connected Bluetooth device.
- If you want to use another Bluetooth device to pair with **goodHearing A1**, switch off the Bluetooth function of the current Bluetooth device, connect the new Bluetooth mobile device to the unit by selecting **goodHearing A1** on the Bluetooth device list and press connect.

- When the Bluetooth device that is connected to the **goodHearing A1** is not playing any sort of media, the device will automatically switch back to being its primary function until the connected device starts playing media or receives a phone call again.
- Should the Bluetooth be disconnected due to range or any other reason, reconnect the **goodHearing A1** with the device that you wish to pair it with.
- If you have connected the **goodHearing A1** to a device and the 2 devices gets out of connecting range, the **goodHearing A1** will have a reminder tone in the background. To switch this tone off, simple turn off the unit for 5 seconds and turn it back on.


5. SAFETY & MAINTENANCE

5.1) Earbud Maintenance

- Clean your ear thoroughly before each use.



- Remove all remaining debris from the earbud and ensure that the hole on the side is also free from any debris after each use.

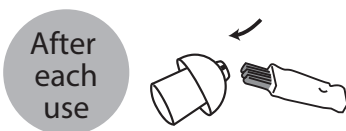


- ※ Use the Cleaning Brush to clean any earwax adhered on the hearing aid, especially in and around the hearing aid and earbud tip.








5.2) Battery Maintenance

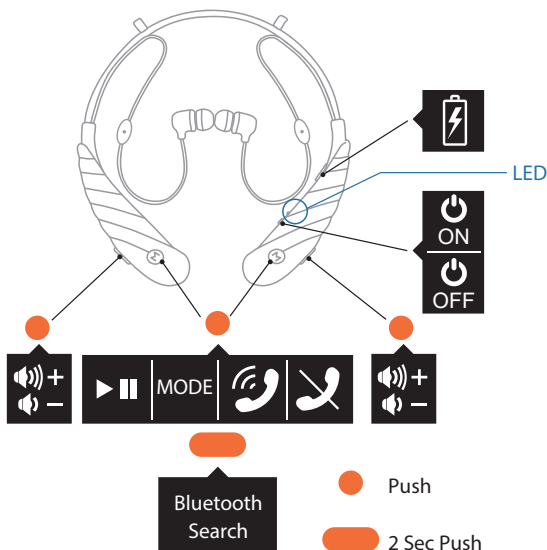- Do not consistently discharge the battery completely before recharging as this may cause depletion of the batteries life.
- We suggest recharging **goodHearing A1** regardless of its battery status before sleeping at night.
- Do not attempt to remove the battery as this may cause damage to the **goodHearing A1** device or may cause injuries.
- Ensure the device is fully charged before removing it from charge.

5.3) LED Indicator status

- Solid RED: Charging
- Solid BLUE: Fully Charged
- Flashing BLUE: ON (working normally)
- Flashing RED and BLUE:
Searching for Bluetooth devices to be paired
- Flashing PURPLE:
Bluetooth Connected with media playing



LED Indicator	
	Red :  ON / OFF
	Blue : Bluetooth Standby
	Red & Blue : Bluetooth Search
	Red-Purple-Blue : Bluetooth Connected with media playing
	Blue : Fully Charged
	Red : Charging

6. FAQ

No Sound ?

- Turn on the device
- Increase the volume
- Clean out the earphone caps
- Check if the earphone cable damages
- Recharge the device

Charging problems?

- Check that the USB is connected correctly
- Ensure that the power supply is normal
(wall outlet switch is on)

Is it compatible for both Android and iOS systems ?

Yes, is compatible for both Android and iOS systems

Can't power on

Out of battery and needs to be charged

There is a lag when I change from listening to media on my phone and back to the hearing aid function is this normal?

Yes, there may be up to 5 second delay depends on the operating system you currently running on your phone

Can I pick up phone calls with the hearing aid ?

Yes, you can, depend on Android or iOS system, your phones will have a default tone that the hearing aid will play as a ringtone

There is still no sound on my SONY phone ?

Please turn off the following sound settings:

Screen lock sound, Vibration on touch,

Keyboard touch tone

There is still no sound on my HTC phone ?

Please turn off the following sound settings:

Touch sounds, Screen lock sounds, Pull-to-

refresh sound, Keyboard touch tones

7. SPECIFICATION

Max. Output @ SPL 90 dB

135 ± 5 dB

Max. Gain

60 ± 5 dB

Continuous Use

Up to 10 hours per charge

Listening Mode

3 Modes

Suitable for

People with moderate to severe hearing loss

Storage Temperature:

14°F ~ 140°F

Storage Humidity

30% ~ 85%

Dimension

160 x 190 x 11 mm

Weight

35 g

Place of Origin

Taiwan

8. WARNING AND REMINDER

- Keep the device dry
- Keep away from direct sunlight
- Keep away from reach of children
- Keep away from high temperature areas
- Avoid liquid near the device
- Avoid discharging the battery for long periods
- Avoid overcharging the device for excessive time periods
- Do not use the device at high volumes
- Do not listen for excessive time periods
- Do not attempt to fix the unit yourself, take it to an authorized dealer
- Ensure the earphones are fitted securely in the ear to reduce sound leak





